

# **Aquafeed Analyzer**

EQUIPMENT SERIES



- Aquafeed and Feed Ingredients
- Moisture, Protein, Fat, Fiber, Starch, Ash and Other Parameters in Under 30 Seconds
- Ready to go Pre-Calibrated with Ingot® Calibrations from Aunir™

# **Process and Quality Control**

Fast, accurate analysis is critical for aquafeed manufacturers to optimize their process and insure quality. Optimizing the process will save the plant time and money and provide very rapid payback. Near Infrared (NIR) analysis is a proven technique designed to provide fast, accurate, and reliable results for the feed and feed ingredients industry.

#### **BENEFITS:**

- Rapid payback
- At-line ready
- · Simple to use
- Flexible sample handling
- Easy data integration
- Low maintenance
- · ISO 12099 Compliant
- Guaranteed calibrations



Spectra Star™

The Unity Aquafeed Analyzer is based on the SpectraStar XT,

a state of the art hi-performance NIR analyzer optimized for the analysis of ground and unground natural products. The

SpectraStar XT contains TAS, a revolutionary technology which monitors the instrument status and precisely calibrates the instrument to first principles. The SpectraStar XT delivers reliable operation and ultimate performance.

The Unity Scientific SpectraStar™ RTW (Rotating Top Window) system gives the ultimate in sampling flexibility and allows the customer to make measurements of all finished feed and feed ingredients.

# Raw material testing to guide proper formulation, to control feed quality, and to optimize least cost formulation

Monitoring and controlling mixers and optimizing the moisture in extruders and dryers

Verifying production batch values, optimizing moisture levels, and documenting finished feed label claims

# Plug and Go Calibrations

The Unity Aquafeed Analyzer is powered by the Ingot® calibration package from Aunir™. This package includes a full suite of calibrations for fresh and salt water aquafeeds, marine meals and other ingredients. Constituents include protein, oil, fiber, ash and starch.

All Ingot calibrations are backed by Aunir's 30 years of experience and are guaranteed to be accurate. Calibration experts are on call to specifically tailor the calibrations to our customer's needs.

POWERED by NGOT Driven by UCal

# **Aquafeed Analyzer**

# Ready to Use

The Unity Aquafeed Analyzer is robust enough to be used both at-line or in the lab. The system features a long life lamp and no air filter or fan, providing for very low cost of ownership. The built in computer is fully networkable and LIMS compatible, allowing for plant management to immediately access results. The Aquafeed Analyzer is a valuable tool to control every step in the feed manufacturing process, from rawingredients to finished feed.

#### The Aquafeed Analyzer package includes the following:

- SpectraStar 2600 XT (1100 nm 2600 nm) with built in 17" touch screen, prediction software, multi-cup adapter, US-RTIR-0016 adapter ring, Unity TAS Standards
- 2 Large cup US-LGOP-0001 for analysis of feed and feed ingredients
- Full Ingot® Calibration package for aquafeeds and ingredients
- One year calibration support and guarantee from Aunir™



#### **HOW TO ORDER**


The Aquafeed Analyzer package includes everything you need to start analyzing finished feeds and ingredients right out of the box.

#### **AQUAFEED ANALYZER**

US-2600-AFD1 SpectraStar 2600 XT (1100 nm - 2600 nm)

#### **AQUAFEED ANALYZER XTR**

US-2600-AFD2 SpectraStar 2600 XT-R (680 nm - 2600 nm)

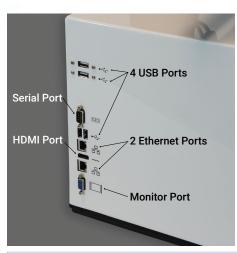

# Specifications Table (Full Specifications in SpectraStar XT Brochure)

- F	
Pre-dispersive scanning monochromator with nominal bandwidth of 10 nm (FWHM)	
1100 - 2600 nm; 680 - 2600 nm XT-R optional	
Tungsten halogen lamp with MTBF rating of 10,000 hours (User changeable)	
Custom dual cooled, extended InGaAs detectors (2)	
10 - 30 seconds	
1.0 nm	
< 0.02 mAu 1100 - 2600 nm	

#### **User Interface**

Built-in Computer	Windows® 7 Embedded, Solid State HDD, 17" high resolution touch screen monitor
Software	UScan™ prediction and data management package
Network Capabilities	OPC Compliant; LIMS Compatible; 4 USB ports, 2 Ethernet (RJ45)

# **Dimensions and Power Requirements**

• • • • • • • • • • • • • • • • • • •	
Size	330 x 381 x 508 mm (with display) 330 x 381 x 270 mm (without display)
Weight	18 kg, 40 lbs.; 22 kg, 48 lbs. with monitor
Voltage	100-240 V AC *, frequency 50-60 Hz, Class 1, protective earth

Unity Scientific Asia Pacific Pty Ltd
Unit 7/ 84-90 Old Bathurst Rd, Emu Plains. NSW, Australia
Ph: 02 4735 1813

E: sydney@unityscientific.com.au www.unityscientific.com.au

Unity Scientific, InfoStar, SpectraStar XL UScan, and UCal are all trademarks of Westco Scientific Instruments, Inc., d/b/a Unity Scientific. All other trademarks are the property of their respective companies. 0901B ©Copyright 2016 All rights reserved. Form #ES\_005\_05.20.16